

## Finecare™ Test Markers

Product Class	Test Item	Diseases	Specimen Type	Reaction Time	Measuring Range	Clinical Reference
Infection	PCT	Sepsis	WB/Serum/Plasma	15 min	0.1-100 ng/mL	0-0.5 ng/mL
	CRP (hsCRP+CRP)	Inflammation	WB/Serum/Plasma	3 min	0.5-200 mg/L	CRP:0-10 mg/L hsCRP:0-1 mg/L
	2 in 1 (CRP+PCT)	Infection & Inflammation	WB/Serum/Plasma	15 min	CRP:0.5-150 mg/L PCT:0.1-100 ng/L	CRP:<10 mg/L PCT:<0.5 ng/L
	SAA	Infection & Inflammation	WB/Serum/Plasma	15 min	1.0-300 mg/L	0-10 mg/L
Coagulation	D-Dimer	Thrombosis	WB/Plasma	5 min	0.1-10 mg/L	0-0.5 mg/L
Cardiac	NT-proBNP	Heart Failure	WB/Serum/Plasma	15 min	18-35000 pg/mL	<75y, 0-300 pg/mL ≥75y, 0-450 pg/mL
	cTn I	Myocardial Infarction	WB/Serum/Plasma	15 min	0.1-50 ng/mL	0-0.3 ng/mL
	Myo	Myocardial Infarction	WB/Serum/Plasma	15 min	2.0-400 ng/mL	0-58 ng/mL
	CK-MB	Myocardial Infarction	WB/Serum/Plasma	15 min	0.3-100 ng/mL	0-5 ng/mL
	H-FABP	Myocardial Infarction	WB/Serum/Plasma	15 min	1-120 ng/mL	0-7 ng/mL
	3 in 1 (cTn I/Myo/CK-MB)	Myocardial Infarction	WB/Serum/Plasma	15 min	Same with single items	Same with single items
	BNP	Heart Failure	WB/Serum/Plasma	15 min	5-5000 pg/mL	0-100 pg/mL
	2 in 1 (cTn I/ NT-proBNP)	Myocardial Infarction	WB/Serum/Plasma	15 min	Same with single items	Same with single items
Tumor	AFP	Liver Cancer	WB/Serum/Plasma	15 min	5-400 ng/mL	0-20 ng/mL
	PSA	Prostate Cancer	WB/Serum/Plasma	15 min	2-100 ng/mL	0-4 ng/mL
	CEA	Cancer	WB/Serum/Plasma	15 min	1-500 ng/mL	0-5 ng/mL
	fPSA	Prostate Cancer	WB/Serum/Plasma	15 min	0.2-30 ng/mL	< 1 ng/mL
Diabetes & Renal Injury	HbA1c	Diabetes Mellitus	WB	5 min	4.0-14.5%	4-6.5%
	MAU	Renal Injury	Urine	3 min	5-300 mg/L	0-20 mg/L
	CysC	Renal Injury	WB/Serum/Plasma	5 min	0.2-10 mg/L	0.5-1.1 mg/L
	NGAL	Renal Injury	Urine	15 min	10-1500 ng/mL	≤131.7 ng/mL
	β2-MG	Proximal Renal Tubule Dysfunction	WB/Serum/Plasma	10 min	0.3-20 mg/L	1.0-2.7 mg/L
Fertility	β-hCG	Fertility	WB/Serum/Plasma	15 min	2-200,000 mIU/mL	5 mIU/mL
	LH	Fertility	WB/Serum/Plasma	15 min	1-100 mIU/mL	Female: 2.95-55.23 mIU/mL Male: 1.70-8.60 mIU/mL
	FSH	Fertility	WB/Serum/Plasma	15 min	1-100 mIU/mL	Female: 4.46-130.00 mIU/mL Male: 1.50-12.40 mIU/mL
	Prog	Fertility	WB/Serum/Plasma	15 min	1.5-60 ng/mL	Female: 7th Week of Pregnancy 24.5±7.6 ng/mL 9-12th Week of Pregnancy 38.0±13.0 ng/mL Male: 0.2-1.5 ng/mL
	PRL	Fertility	WB/Serum/Plasma	15 min	1-200 ng/mL	Female: 4.60-25.07 ng/mL Male: 3.45-17.42 ng/mL
Hormone	T3	Thyroid Malfunction	WB/Serum/Plasma	15 min	0.61-9.22 nmol/L	1.3-3.1 nmol/L
	T4	Thyroid Malfunction	WB/Serum/Plasma	15 min	12.87-300 nmol/L	66-181 nmol/L
	TSH	Thyroid Malfunction	WB/Serum/Plasma	15 min	0.1-100 mIU/L	0.3-5.2 mIU/L
	Vitamin D	Bone Metabolism	Plasma/Serum	15 min	5-100 ng/mL	30-100 ng/mL

# Wondfo®

## Finecare™ FIA Meter Plus (FS-113)

Fluorescence Immunoassay Rapid Quantitative Test

- Rapid
- Accurate
- Convenient
- Reliable
- Popular
- Extensible



Wondfo®

Guangzhou Wondfo Biotech Co., Ltd.  
 Add: No. 8 Lizhishan Road, Science City, Luogang District, 510663, Guangzhou, P.R. China  
 Tel: +86-20-32296083 Fax: +86-20-32296063  
 E-mail: sales@wondfo.com.cn Websites: www.wondfo.com.cn

For Medical Supplies  
 01010060298 - 01210629090  
 www.delta-care.com

دلتا كير  
 DELTA CARE

الوكيل الوحيد  
 بجمهورية مصر العربية  
 شركة

## Fineware™ FIA Meter Plus

Fineware™ FIA Meter Plus is a fluorescence immunochromatographic analysing system which can help diagnose conditions such as infection, diabetes, cardiovascular diseases, renal injury and cancers, etc.

The Fineware™ platform provides reliable and quantitative results of various kinds of analytes in human blood or urine within several minutes.

**Application:** Laboratory, ER, Cardiology, ICU, Respiratory, Pediatrics, etc.

### Test Item

#### ● Fertility

β-HCG, LH, FSH, PRL, Prog, E2, Testosterone

#### ● Infection

PCT, CRP (hsCRP+CRP), 2 in 1 (CRP+PCT), SAA

#### ● Cardiac

NT-proBNP, cTn I, Myo, CK-MB, H-FABP, BNP  
3 in 1 (cTn I + Myo + CK-MB),  
2 in 1 (cTn I+NT-proBNP)

#### ● Diabetes & Renal Injury

HbA1c, MAU, Cys C, NGAL, β2-MG

#### ● Tumor

AFP, PSA, CEA, fPSA

#### ● Coagulation

D-Dimer

#### ● Thyroid Hormone

T3, T4, TSH, FT4, FT3


#### ● Hormone

Vitamin-D  
Cortisol

#### ● Ongoing

Ferritin,  
S100β,  
ST2,  
Vitamin B12,  
H. pylori,  
AMH

### Principle and Operation

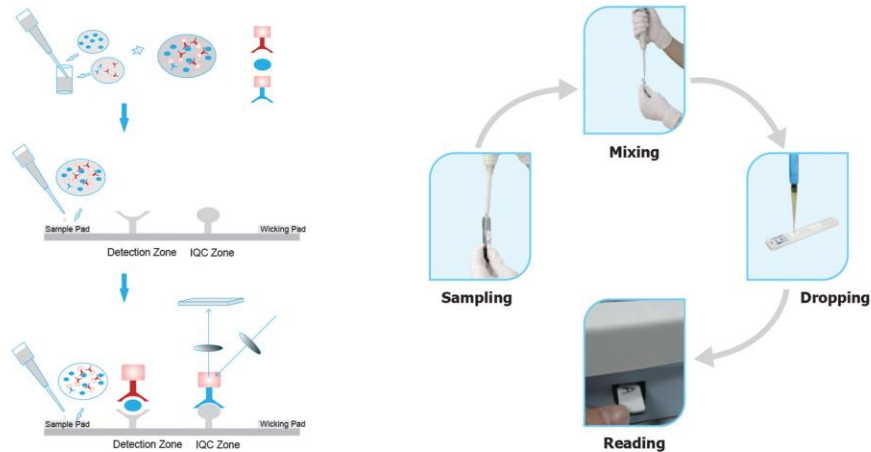

## Why Choose Fineware™



Extra quick operating mode for large number of specimens

Unique quick mode



Specially-designed sampler for small volume of specimen


Designed sampler



Powerful information management (Direct connection to LIS/HIS)

Information management

Reliable accuracy

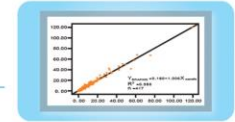

Good correlation with wet chemistry method results

Instant result

Room temperature storage



Solid synthetic craft, can be stored at room temperature (4-30°C)



Convenient and fast. Get results in 3-15 minutes

دلتا کیر  
DELTA CARE

Built-in Thermal Printer

Real-time printing

8" LCD Touch Screen

Inner Temperature Control



**Dimensions:**  
270(L) x 238(W) x 146(H)mm  
**Weight:**  
2.8kg  
**Power adapter:**  
100-240AC, 50-60Hz

LIS Connection  
Direct connection to LIS / HIS

ID Chip Port  
The ID chip contains all information of test item and lot number

Test Channel  
The Test Cartridge Holder